VULPES GOODS
BODYCARE

RED LIGHT THERAPY ELITE

MANUAL - HANDLEIDING

WWW.VULPESGOODS.COM

ENGLISH

4

The Device

3

User Guide

7

Characteristics and Application

9

Product parameters

11

NEDERLANDS

12

Het apparaat

13

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16

Kenmerken en Toepassingen



18

Product parameters

20

ENGLISH

THE DEVICE

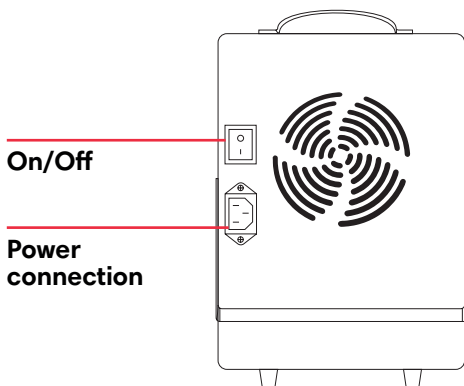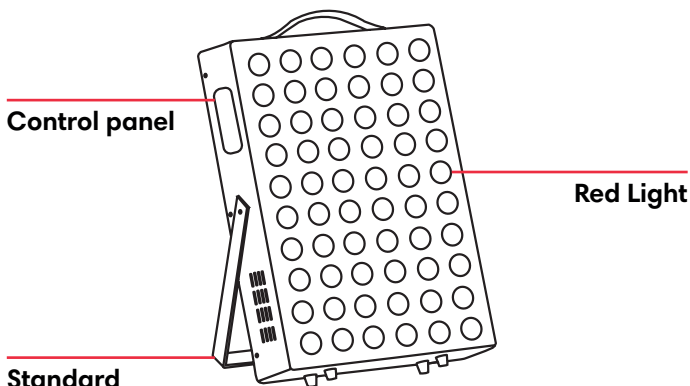Red Light

Position the red lights at the desired location before switching on the device. Only turn on the lights once everything is correctly installed. Never look directly into the lights during use. Use the supplied sunglasses during use.

Power Connection

Connect the power cable to the power connection to supply the device with power. Without the power connection, the Red Light will not operate.

On/Off

Set the switch to I to turn the Red Light on. Set the switch to O to turn the Red Light off.

Standard


Use the stand to position the Red Light at the desired angle for use.

Control Panel

Use the control panel to turn on the Red Light, switch between modes, adjust the brightness, or set a timer. Please refer to the next page for a detailed explanation of the control panel.

Remote Control

The Red Light can be operated using the included remote control. The buttons on the remote control correspond to the buttons on the Red Light. It is not possible to activate the 10-minute timer using the remote control.

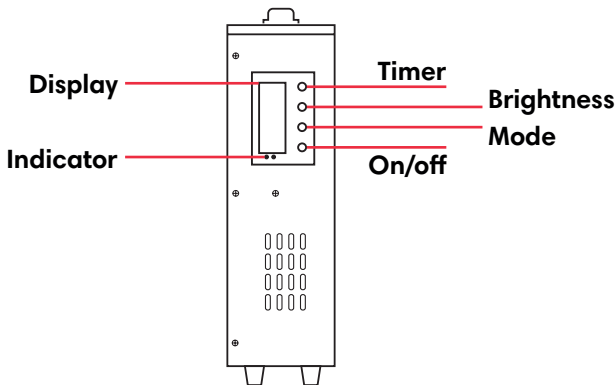Timer

Press the timer button once to activate a 10-minute timer. Each additional press adds 10 minutes, up to a maximum of 60 minutes.

Brightness

Press the brightness button once to increase the brightness. The brightness can be adjusted up to level 5.

Mode

The Red Light features three modes: 850 nm near infrared (invisible), 650 nm red light (visible), and 850 nm + 660 nm (Dual Light Mode). Each press of the mode button switches to the next mode. The selected mode is shown on the indicator.

Indicator

Use the indicator to identify the selected mode via the illuminated lights below the display.

Display

Use the display to view the selected mode, brightness level, and remaining time.

USER GUIDE

1. Place the Red Light on a stable surface and connect the power cable. Plug the adapter into a power outlet.
2. Set the power switch to I and press the On/Off button to turn on the device. The device automatically starts in the default 850 nm (NIR) mode.
3. Adjust the angle of the lamp so that the light is directed at the desired treatment area.
4. Position the treatment area at a distance of approximately 10–15 cm from the light source. Keep sufficient distance to avoid discomfort or burns.
5. Use the control panel to select the desired mode, brightness level, and treatment time.
6. Start the treatment session. The recommended duration is 10 to a maximum of 30 minutes per session.
7. The device will automatically switch off when the set time has elapsed. Do not exceed the maximum usage time of 30 minutes.
9. After use, turn off the device and unplug the power adapter.

Mode

The Red Light features multiple light modes, each designed for specific applications. Select the mode that best matches your intended use.

Near Infrared – 850 nm (invisible light)

This mode is designed for deeper tissue exposure and general wellness use. Near infrared may support relaxation and recovery and is commonly used for full-body applications.

Red Light – 660 nm (visible red light)

This mode is suitable for surface-level use, such as skin care. Red light is often used to support skin condition and cosmetic applications.

Dual Light Mode – 850 nm + 660 nm

The dual light mode combines near infrared and red light. This mode is suitable when both surface-level and deeper exposure are desired and is commonly used for combined applications.

Always select the appropriate mode based on personal comfort and follow the recommended usage time.

CHARACTERISTICS AND APPLICATION

- Avoid direct exposure of the eyes to red light. Use eye protection such as goggles or an eye mask.
- Stop using the device immediately if redness, itching, pain, or discomfort occurs and consult a medical professional.
- Do not move the lamp or touch the device during operation to prevent accidents.
- This device is not suitable for children under 13 years of age or for pregnant women.
- Do not use the device while sleeping, unconscious, or under the influence of alcohol.
- Use with caution if you have sensitive skin or sensitivity to sunlight or red light.
- Keep the device and accessories away from water, heat, and fire. Do not immerse the device in water or other liquids.
- Do not touch the device or power cord with wet hands.
- Do not operate the device if the power cord, plug, or housing is damaged.
- After use, always unplug the device. Never leave it unattended while plugged in to avoid electrical hazards.
- Do not use the device on open wounds, bleeding areas, or damaged skin.



- Do not exceed a maximum usage time of 30 minutes per session.
- Do not attempt to disassemble the device. Unauthorized modification or repair will void the warranty.
- Remove all metal jewelry and accessories from the treatment area before use. Do not use over clothing with metal parts such as zippers or buttons.
- After each use, turn off the device and unplug the adapter to prevent damage or safety risks.
- Regularly clean the lamp surface with a soft, dry cloth to maintain optimal performance.
- Protect the device from impacts, drops, and pressure to prevent damage.

PRODUCT PARAMETERS

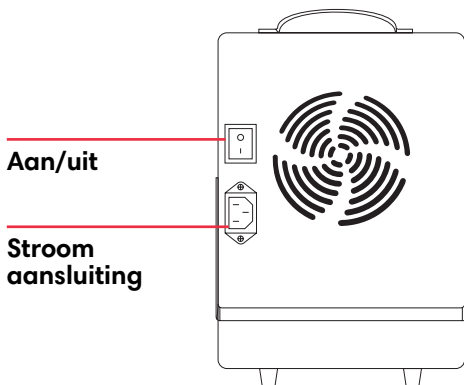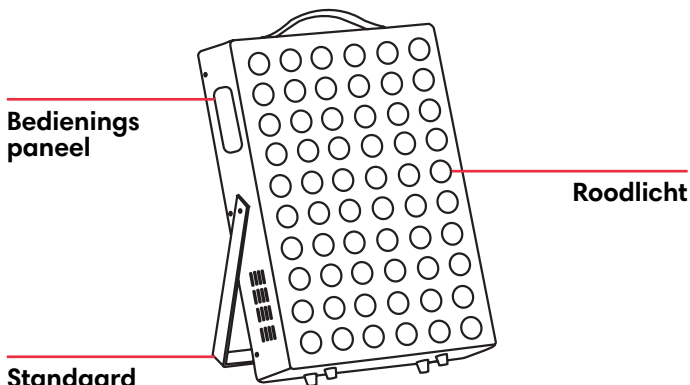LED Power:	150 W
Material:	Aluminum housing, PVC, and stainless steel
LED Quantity:	60 pcs
Input Voltage:	AC 100–245 V / 50–60 Hz
Product Size:	270 × 192 × 65 mm
Timer Settings:	20 / 30 / 40 / 50 / 60 minutes
Operating Temperature:	-20 °C to 40 °C
Waterproof Rating Recommended:	IP3
Irradiation Distance:	10–30 cm
Dimming Levels:	Level 1–5
Spectral Band:	660 nm + 850 nm

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