VULPES GOODS
BODYCARE

RED LIGHT THERAPY PRO V2

MANUAL - HANDLEIDING

WWW.VULPESGOODS.COM

ENGLISH

3

The Device	4
User Guide	6
Characteristics and Application	9
Product parameters	10

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
THE DEVICE

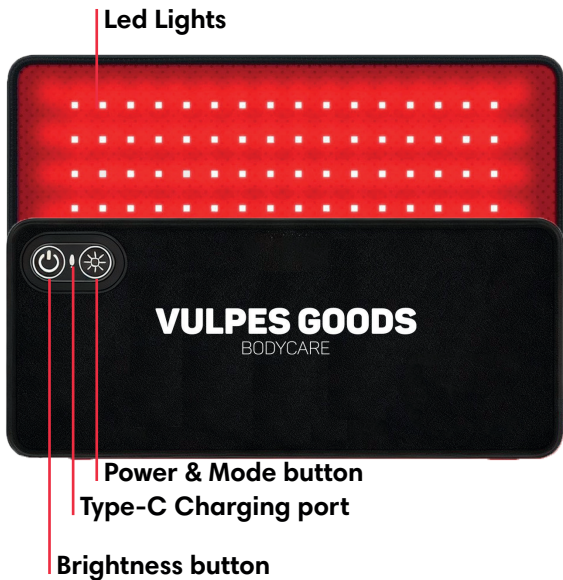LED Lights

The LED lights emit red light and near-infrared light, depending on the selected mode. Place the LED side on or above the desired treatment area and use the device only on a dry, flat surface.

Power & Mode Button

Press and hold the button to turn on the device. The device starts in Normal Mode by default. Short-press the button to switch between the available modes. Press and hold the button again to turn off the device.

Type-C Charging Port

Use the Type-C charging port to charge the device. During charging, the red indicator light flashes. When the battery is fully charged, the indicator light turns green. It is recommended to use the device when the battery is fully charged.

Brightness Button

Press the brightness button to adjust the light intensity. Press repeatedly to cycle through the available brightness levels and select a level that feels comfortable.

USER GUIDE

1. Make sure the battery is fully charged before using the device. During charging, the red indicator light flashes. When charging is complete, the light turns green and remains steady.

2. Press and hold the power button to turn on the device. The device starts in Normal Mode by default.

3. Short-press the power button to switch between modes.

First short press: Flicker Mode 1

Second short press: Flicker Mode 2

The modes cycle in this order.

4. Press the brightness button to adjust the light intensity. Repeat to switch between the available brightness levels.

5. Press and hold the power button again to turn off the device.

6. For safety reasons, the device automatically shuts off after 30 minutes of use.

Red + NIR Mode (660 nm + 850 nm)

In this mode, the device emits both 660 nm red light and 850 nm near-infrared light. This combined output provides dual light exposure and is intended for general use. This mode is commonly selected to support relaxation, ease tension, and promote overall wellness.

NIR Mode (850 nm)

In this mode, the device emits only 850 nm near-infrared light. This mode is suitable for targeted use and is often chosen to support a calming effect and localized comfort.

Red Mode (660 nm)

In this mode, the device emits 660 nm red light. This mode is suitable for surface-level use and is commonly used to support a gentle, revitalizing sensation and help maintain the body's natural balance.

Skin Test and Usage Advice

Before first use, it is recommended to perform a simple skin sensitivity test to observe how your skin reacts to the light.

1. Place the device on your arm, switch it on, and use it for 5 to 15 minutes. If your skin becomes noticeably red or experiences strong itching, this may indicate light sensitivity. In this case, stop using the device and consider consulting a specialist before further use.
2. Before use, clean the treatment area and ensure the skin is free from makeup or cosmetic products. Some ingredients may increase sensitivity to light.
3. After use, mild redness or slight itching may occur. This is usually temporary and should subside within 1 to 2 hours. If discomfort persists, a burning sensation occurs, or other symptoms develop, stop using the device immediately.

Note:

Perform the skin test and use the device under normal lighting conditions. Discontinue use immediately if you feel unwell or experience sensitivity to the light.

CHARACTERISTICS AND APPLICATION

- The Red Light Pad is designed for controlled and comfortable body use.
- The product features multiple light modes and an integrated temperature protection system to support safe operation.
- The pad is flexible but must always be used flat and must not be folded.
- The device is not water-resistant and is intended for use in dry environments only.
- Use the Red Light Pad on a stable, dry surface and keep it away from water and moisture.
- Do not use the product near heat sources or heating appliances, such as electric blankets.
- If the LED indicators flash rapidly, the temperature protection is activated and use should be paused until the product has cooled down.
- The Red Light Pad is suitable for general wellness and relaxation purposes.
- The product is intended for use on intact skin and only while the user is awake and alert.
- Use the Red Light Pad indoors and always follow the instructions provided in the user manual.

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