VULPES GOODS
BODYCARE

RED LIGHT THERAPY STANDING

MANUAL - HANDLEIDING

WWW.VULPESGOODS.COM

ENGLISH

4

The Device

3

User Guide

7

Characteristics and Application

8

Product parameters

10

NEDERLANDS

11

Het apparaat

12

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15

Kenmerken en Toepassingen


16

Product parameters

18

ENGLISH

THE DEVICE

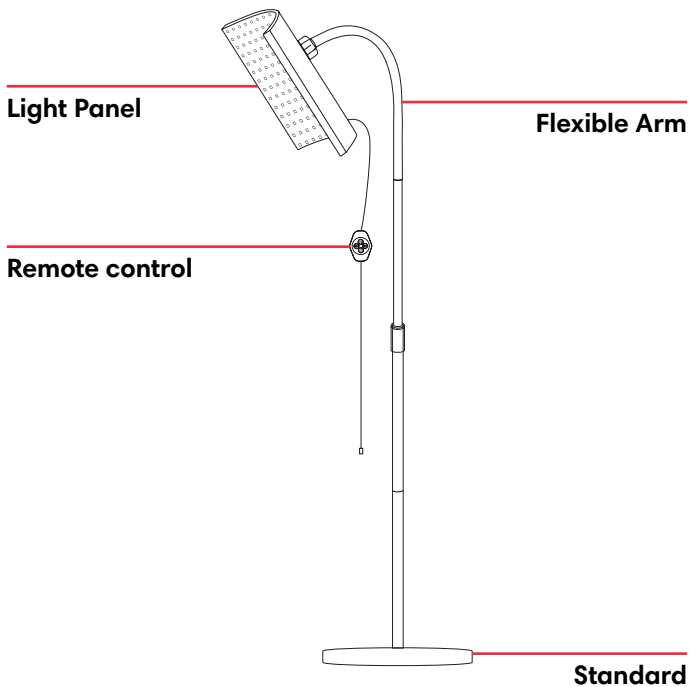Light Panel

The light panel contains the LED lights and is attached to the flexible arm. Aim the panel at the desired treatment area and ensure it is securely positioned during use.

Flexible Arm

The flexible arm allows you to easily position the light panel. Gently bend the arm to direct the light accurately.

Remote Control

The remote control allows you to operate the device, including turning it on or off, adjusting brightness and setting the timer.

The separate remote control allows full operation of the device. The device can be turned on or off, brightness can be adjusted, a mode can be selected, and a timer can be set.


Use the ON/OFF buttons to power the device on or off. The brightness buttons are used to increase or decrease the light intensity.

The mode button (M) switches between constant light, low-frequency flicker mode, and high-frequency flicker mode.

The timer buttons set an automatic shut-off of 20, 30, 40, 50, or 60 minutes. Press once to start the timer and press again to cancel it. The selected time is confirmed by a light signal.


Standard

The stand provides stable placement of the device on the floor. Please see the next page for proper installation and assembly of the stand.

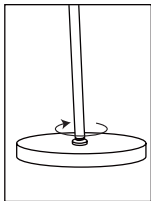
Step 1: Attach the pillar to the base

Insert the pillar into the base and tighten it until the pillar stands firmly and upright.

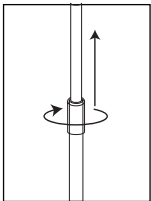
Step 2: Adjust the height

Pull up the telescopic rod to the desired height and tighten the locking ring to secure the position.

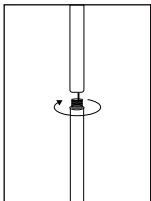Step 3: Connect the flexible arm to the pillar

Screw the flexible arm securely onto the top of the pillar to ensure stability.

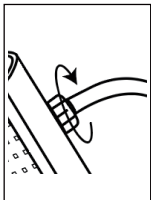Step 4: Attach the light panel

Attach the light panel to the flexible arm and tighten it. Adjust the panel to the desired position for use.

USER GUIDE

1. Turn on the device using the power switch. The lamp will begin to warm up; wait 1–2 minutes until the light becomes stable.
2. Adjust the height and angle of the lamp so the red light is directed at the desired treatment area.
3. Maintain a distance of approximately 10–15 cm between the lamp and the skin to avoid overheating or discomfort.
4. Start the treatment session. Each session should last between 10 and a maximum of 30 minutes. Do not exceed this duration.

Usage Frequency

1. Use the lamp 2–3 times per week.
2. After a period of consistent use, the frequency may be adjusted based on personal comfort and experience.

CHARACTERISTICS AND APPLICATION

- Avoid direct exposure of the eyes to red light. Use eye protection such as goggles or an eye mask.
- Stop using the device immediately if redness, itching, pain, or discomfort occurs and consult a medical professional.
- Do not move or touch the device during operation to prevent accidents.
- The device is not suitable for children under 13 years of age or for pregnant women.
- Do not use the device while sleeping, unconscious, or under the influence of alcohol.
- Use with caution if you are sensitive to sunlight or red light.
- Keep the device and accessories away from water, heat, and fire. Do not immerse the device in any liquids.
- Do not touch the device or power cord with wet hands.
- Do not use the device if the power cord, plug, or housing is damaged.
- After use, always unplug the device. Never leave it unattended while plugged in.
- Do not use the device on open wounds, bleeding areas, or damaged skin.


- Do not exceed a maximum usage time of 30 minutes per session.
- Do not attempt to disassemble or modify the device. Unauthorized repairs will void the warranty.
- Remove all metal jewelry and accessories from the treatment area before use. Do not use over clothing with metal parts.
- After each use, turn off the device and unplug the adapter to prevent damage or safety risks.
- Clean the surface of the device regularly with a soft, dry cloth.
- Protect the device from impacts, drops, and pressure to avoid damage.

PRODUCT PARAMETERS

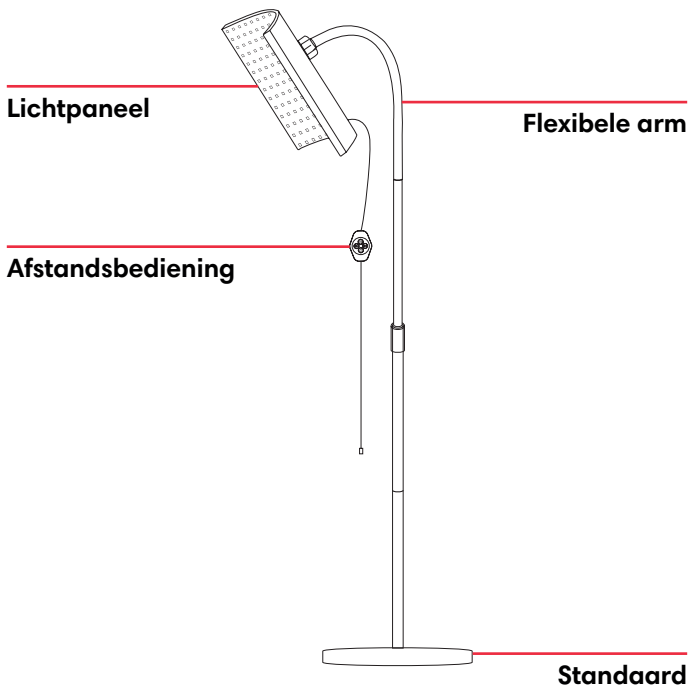Power:	21 W
Material:	Aluminum housing, PVC, and stainless steel
Installation Method:	Base / telescopic pole / light panel
Power Connector:	DC5521 interface
Adapter:	DC 12 V / 2 A
Timer Settings:	20 / 30 / 40 / 50 / 60 minutes
Adjustable Height:	35–150 cm
Operating Temperature:	–20 °C to 40 °C
Waterproof Rating:	IP45
Recommended Irradiation Distance:	10–30 cm
Dimming Levels:	Level 1–5
Spectral Bands:	660 nm + 850 nm

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