VULPES GOODS
BODYCARE

RED LIGHT THERAPY MAX

MANUAL - HANDLEIDING

WWW.VULPESGOODS.COM

ENGLISH

3

The Device

4

User Guide

7

Characteristics and Application

9

Product parameters

11

NEDERLANDS

12

Het apparaat

13

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16

Kenmerken en Toepassingen


19

Product parameters

21

ENGLISH

THE DEVICE

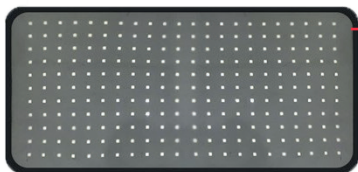LED-lights

Power port

Strap

Power adapter

Controller

LED Lights

The LED lights are integrated into the mat and provide light output during use. Place the mat flat on the desired body area so the LED lights maintain good contact for comfortable use.

Power Port

The power port is used to connect the power adapter. Make sure the adapter is securely connected before turning on the device.

Strap

The adjustable strap is used to keep the mat in place during use. Fasten the strap firmly but comfortably to prevent the mat from shifting.

Power Adapter

The power adapter supplies electricity to the device. Use only the supplied adapter and ensure it is properly connected to a power outlet.

Controller

With the controller, the device can be operated. The controller is used to turn the device on and off and to adjust the mode, power level, and time settings

Mode Selection

Use the mode button to switch between Regular Mode, 10 Hz pulse mode, and 40 Hz pulse mode. Each mode is selected by short-pressing the button.

Power Adjustment

Use the plus and minus buttons to adjust the power level from P1 (low) to P5 (high). Always select a level that feels comfortable.

Time Setting

Use the time buttons to set the session duration between 10 and 90 minutes. The device switches off automatically when the set time ends.

USER GUIDE

Using the red light therapy mat is simple and suitable for daily use. Each session should last up to 30 minutes per use, depending on personal comfort.

1. Connect the power cord to the red light therapy mat and plug the adapter into a power outlet. Make sure all connections are properly secured before turning on the device.

2. Press and hold the power button on the controller for approximately 2 seconds to turn on the device. The device starts by default in Regular Mode (P1) with a preset time of 30 minutes.

3. After powering on, the settings can be adjusted as needed:

Mode Selection (ON/OFF button):

Short press to cycle through:

Regular mode

10 Hz pulse mode

40 Hz pulse mode

Power Level (PWM + / -):

Use the plus and minus buttons to adjust the intensity level from low (P1) to high (P5). Select a level that feels comfortable.

Time Setting (TIME + / -):

Set the session duration between 10 and 90 minutes. Increase or decrease the time gradually using the plus or minus buttons.

4. Place the red light therapy mat on the desired body area and start the session. A duration of 20–30 minutes is generally recommended for comfortable use. For larger or less sensitive areas, longer sessions may be selected, provided comfort is maintained.

5. When the set time ends, the device switches off automatically. Unplug the device after use and store it in a cool, dry place.

CHARACTERISTICS AND APPLICATION

- Store the red light therapy mat in a cool and dry place when not in use.
- Do not allow children to use the red light therapy mat without adult supervision.
- Do not use this device on infants or young children.
- Do not allow sleeping, unconscious, or non-self-reliant persons to use the device.
- Do not exceed a maximum usage time of 90 minutes per session.
- Do not use the device on open, fresh, or damaged wounds.
- Do not use the device together with creams, ointments, or products containing heat-stimulating or irritating ingredients.
- Do not disassemble, modify, or attempt to repair the device. There are no user-serviceable parts inside.
- Do not use the device in or near water or other liquids.
- Do not drop the device into water or place it where it could fall into liquids.
- Do not touch the device, plug, or switches with wet hands.
- Do not apply or touch the device to wet or damp skin.


- Do not pull, carry, or lift the device by the power cord.
- Do not use the device on or near heated surfaces.
- Do not operate the device in areas where flammable or combustible substances or vapors may be present.
- Do not store the device at temperatures exceeding 50 °C.
- Do not use the device if it is damaged. Continued use of a damaged device may result in injury or electric shock.
- Do not use extension cords unless they are rated to carry the total power load of the device.
- Always unplug the device by grasping the plug, not the cord.
- Unplug the device when not in use to avoid tripping hazards.

PRODUCT PARAMETERS

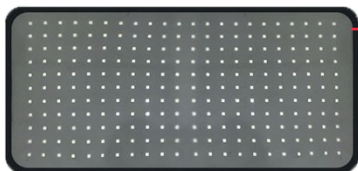Lamp Power:	minimum 10 W, maximum 35 W
Irradiance:	85–114 mW/cm ²
Wavelengths:	660 nm and 850 nm
Input Voltage:	110–240 V
Dimensions:	approx. 61 × 30.5 cm
Material:	Neoprene fabric and leather
Operation Modes:	Regular mode 10 Hz pulse mode 40 Hz pulse mode

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